

Laryngeal FS Rhodium Mirrors




Laryngeal FS Rhodium mirrors are coated with Rhodium on the front surface. Rhodium is a precious metal as gold and as a precious metal it is chemical resistant, even against acids.

The features are:


Images of exact sharpness and contrast up to the smallest detail, without any distortion

only one reflection
eliminates the so-called »ghost image«
→ exact sharp image

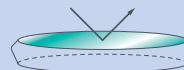

Front Surface Mirror

Two reflections
→ double image

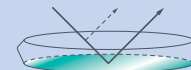


Rear Surface Mirror

Bright and brilliant reflection in natural colours
thanks to an optimal reflection spectrum.

Highest sight comfort
gentle to the eye and good for fatigue-proof working

Easy cleaning
due to the anti-stick effect of the special
HAHNENKRATT front surface coating

Unlimited resistance
against disinfection and sterilisation, even against
acid containing refresh liquids.

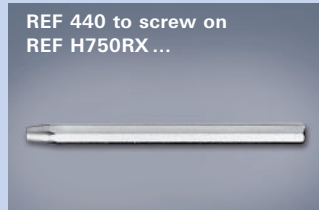

Laryngeal Mirrors Rhodium plane	Ø in mm	REF without thread	REF with thread	REF Mirror with handle hexagonal	REF Mirror with handle round
each	6	H951RX06	H950RX06	H956RX06	H957RX06
	8	H951RX08	H950RX08	H956RX08	H957RX08
	10	H951RX10	H950RX10	H956RX10	H957RX10
	12	H951RX12	H950RX12	H956RX12	H957RX12
	14	H951RX14	H950RX14	H956RX14	H957RX14
	16	H751RX16	H750RX16	H756RX16	H757RX16
	18	H751RX18	H750RX18	H756RX18	H757RX18
	20	H751RX20	H750RX20	H756RX20	H757RX20
	22	H751RX22	H750RX22	H756RX22	H757RX22
	24	H751RX24	H750RX24	H756RX24	H757RX24
	26	H751RX26	H750RX26	H756RX26	H757RX26
	28	H751RX28	H750RX28	H756RX28	H757RX28
	30	H751RX30	H750RX30	H756RX30	H757RX30

Laryngeal Mirrors Handles		REF	REF		
each		400	440		