

MEGA FS

FRONT SURFACE


The MEGA FS is the alternative to the MEGA FS Rhodium. Rhodium is even more expensive than gold. In order to offer a more favourable front surface mouth mirror, we developed a unique combination of reflection layers on the base of our decades-long experience.

Both mouth mirrors distinguish in the way of reflection. Rhodium reflects colours with a »warmer« touch, the new FS-Combi coating reflects with a »cooler« touch.




Images of exact sharpness and contrast up to the smallest detail, without any distortion

only one reflection
→ exact sharp image

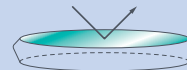

Front Surface Mirror

Two reflections
→ double image

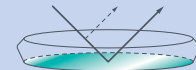

Rear Surface Mirror



Bright and brilliant reflection in natural colours
thanks to an optimal reflection spectrum.

High sight comfort

gentle to the eye and good for fatigue-proof working.

Easy cleaning

due to the anti-stick effect of the special HAHNENKRATT front surface coating

Unlimited resistance

against disinfection and sterilisation, even against acid containing refresh liquids.



| MEGA FS mouth mirrors | | REF | |
|-------------------------|---------|-----------|-----------|
| Box of 12 size | Ø in mm | SS-thread | CS-thread |
| | | plane | plane |
| 4 | 22 | 6 560 | 6 565 |
| 5 | 24 | 6 561 | 6 566 |