

# MEGA FS FS

## FRONT SURFACE


The MEGA FS is the alternative to the MEGA FS Rhodium. Rhodium is even more expensive than gold. In order to offer a more favourable front surface mouth mirror, we developed a unique combination of reflection layers on the base of our decades-long experience.

Both mouth mirrors distinguish only in the way of reflection. Rhodium reflects colours with a »warmer« touch, the new FS-Combi coating reflects with a »cooler« touch.




Images of exact sharpness and contrast up to the smallest detail, without any distortion

only one reflection =  
exact sharp image

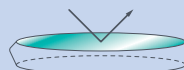

Front Surface Mirror

Two reflections =  
double image

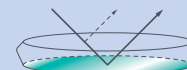

Rear Surface Mirror



**Bright and brilliant reflection in natural colours**  
thanks to an optimal reflection spectrum.

**High sight comfort**  
gentle to the eye and good for fatigue-proof working.

**Easy cleaning**  
due to the anti-stick effect of the special  
HAHNENKRATT front surface coating

**Unlimited resistance**  
against disinfection and sterilisation, even against  
acid containing refresh liquids.



MEGA FS   mouth mirrors		REF	
Box of 12 size	Ø in mm	SS-thread	CS-thread
		plane	plane
4	22	6 560	6 565
5	24	6 561	6 566