

TOPvision FS

FRONT SURFACE


The TOPvision FS is the alternative to the TOPvision FS Rhodium. Rhodium is even more expensive than gold. In order to offer a more favourable front surface mouth mirror, we developed a unique combination of reflection layers on the base of our decades-long experience.

Both mouth mirrors distinguish only in the way of reflection. Rhodium reflects colours with a »warmer« touch, the new FS-Combi coating with a »cooler« touch.




Images of exact sharpness and contrast up to the smallest detail, without any distortion

only one reflection =
exact sharp image

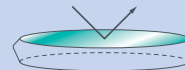

Front Surface Mirror

Two reflections =
double image

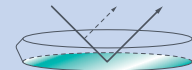

Rear Surface Mirror



Bright and brilliant reflection in natural colours
thanks to an optimal reflection spectrum.

High sight comfort


gentle to the eye and good for fatigue-proof working.

Easy cleaning

due to the anti-stick effect of the special
HAHNENKRATT front surface coating

Unlimited resistance

against disinfection and sterilisation, even against
acid containing refresh liquids.

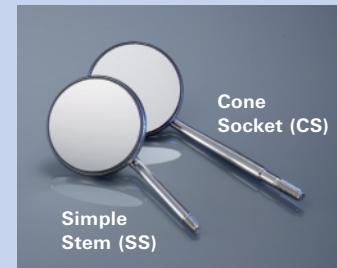

TOPvision FS mouth mirrors		REF	
Box of 12 size	Ø in mm	Simple Stem (SS)	Cone Socket (CS)
		plane	plane
4	22	633X4	636X4
5	24	633X5	636X5