

SEplus

MOUTH MIRRORS

HIGH QUALITY IN EACH DETAIL

Since decades dentists all over the world rely on the superior quality of our SEplus mouth mirrors.

In the dental practice mouth mirrors are daily sterilised and disinfected several times. Considering these extreme conditions, the superior quality – in every single detail – is decisive.



Best reflection and image quality


The SEplus reflection layer with its fivefold pure silver coating guarantees an extreme high reflection. The result is a bright and most brilliant reflection in natural colours. Thus the SEplus provides a high sight comfort.

Best stability against corrosion

According to our experience of over 50 years, the stainless steel is stable against corrosion. It complies with ISO 7153-1 (Surgical Instruments – Metallic materials).

Best resistance against disinfection and sterilisation

The silver coating is protected by a combination of layers, mainly made out of copper and titanium, which was developed by our company. Additionally the mirror is sealed with a special protective film. This unique combination ensures the extreme long life of SE plus mouth mirrors.

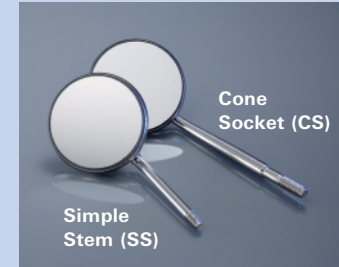

SEplus mouth mirrors		REF			
Box of 12 size	Ø in mm	Simple Stem (SS)		Cone Socket (CS)	
		plane	magnifying	plane	magnifying
3	20	702X3	706X3	711X3	716X3
4	22	702X4	706X4	711X4	716X4
5	24	702X5	706X5	711X5	716X5